




For rapid sterilization

**EQUITRON®**

## TABLE TOP STERILIZER with drying

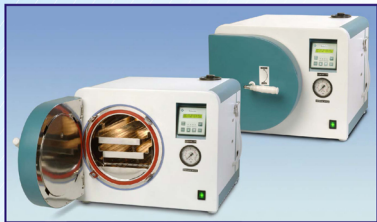

### Applications:

- ♦ EQUITRON Table Top High Speed Sterilizers are ideal for all applications requiring rapid, reliable and total destruction of all living microorganisms.
- ♦ Ideal for Ophthalmic Surgeons, Dental Surgeons, Surgical Doctors, Polyclinics and Hospitals.

### EQUITRON Advantages:

- ♦ Sterilizing, exhaust of water and steam drying are fully automatic.
- ♦ Requires a sterilizing time of just 7 minutes at 134°C (instead of 20 minutes at 121°C).
- ♦ A totally hassle and tension free sterilizer which can even be kept inside an operation theatre.
- ♦ Results reproducible because full cycle is automatic, requiring no human intervention.

### Construction:

- ♦ Full stainless steel construction.
- ♦ Single lever pronged closure system with patented Silicone Gasket (registered design No. 208774.)

### Standard Features:

- ♦ µP Controller with following parameters programmable by user: **Sterilizer Time:** up to 99 minutes. **Sterilizing Temperature:** up to 134°C. **Drying time:** up to 99 minutes.
- ♦ Pt100 Sensor for precise monitoring.
- ♦ End cycle buzzer and auto reset.
- ♦ Valves for water filling and exhaust.
- ♦ Industrial Grade energy efficient heater.
- ♦ Drying with door open facility.

### Safety:

- ♦ Sensor open alarm and over temperature protection.
- ♦ Pressure interlock on door.
- ♦ Spring loaded safety valve for over pressure.

| Model    | Chamber Size            | External Dimension<br>W/D/H in mm | Tray    | Heater | App. Shipping Details (Air)<br>Net / Gr / Volume |
|----------|-------------------------|-----------------------------------|---------|--------|--------------------------------------------------|
| #7471-16 | Ø 230 x 400 mm (16 lit) | 560 x 500+125 x 450               | 2 trays | 1.2 kW | 50 kg / 65 kg / 0.3 m <sup>3</sup>               |

**MultiLab**

No.51/23, Arathoon Road, Off Appavu Lane,  
Royapuram, Chennai - 600 013. Tamil Nadu, India  
Ph:044-25984175, 45526026, 094441 62346  
Email: multilabindia@gmail.com, www.multilab.biz