

Shear/Scratch Tester 550/551

TABER
INDUSTRIES
ISO 9001:2000 CERTIFIED

Surface Finish Testing

Precision instrument measures damage to surfaces not classified as ordinary wear.


Using interchangeable tools, the Taber Shear / Scratch Tester is used to measure the susceptibility of material surfaces to scratching, marring, shearing, gouging, scraping, engraving, and other similar physical damage not classified as ordinary wear.

Established tester can be used in diverse applications.

Useful in quality control as well as materials and product development, the applications are diverse. The Shear / Scratch Tester is ideal for evaluating adhesive qualities of surface coatings, detecting the lack of homogeneity in comparative materials, rating similar or competitive materials, or analyzing the effects of age under controlled environmental conditions.

Test a wide range of products.

Precision, calibrated weights provide a flexible testing range for a variety of products. The Shear / Scratch Tester can test rigid organic materials, adhesive coatings, paints, powder coatings, anodized materials, soft metals, synthetics, linoleum, laminates, plastics, glass, plus many others.

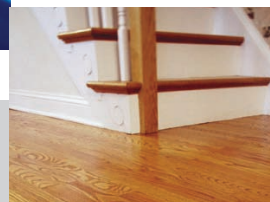

Shear/Scratch Tester 550/551

Surface Finish Testing

Available In Two Models:

Model 550 - 115V, 60Hz

- Turntable speed 0.5 rpm
- Conical diamond tool with 90°, 3 mil radius point

Model 551 - 230V, 50Hz 

- Turntable speed 5 rpm
- Conical diamond tool with 90°, 3.5 mil radius point

Standard Features:


Constant speed motor
Balanced, calibrated scale beam
Variable load from 0 - 1000 grams
Adjustable gage block tests up to 1/2" (12.5mm) thick specimens
Soft touch ON/OFF switch
Interchangeable tools

Optional Accessories (sold separately):

Contour Shear Tool (Honda and Toyota requirements)
Diamond Scratch Tool (converging point 3 sides of cube)
Scratch Tip, 1.0 mm diameter hemisphere
Diamond Standardization Plates
Scale Beam Height Extension Kit
15N Scale Beam Fixed Weight

Shear / Scratch Tester Includes:

500 Gram Calibrated Weight (2 each)
S-20 Precision Ground Tungsten Carbide Tool
Conical Diamond Tool
10x Measuring Magnifier
Spirit Bubble Level
Hex Wrench

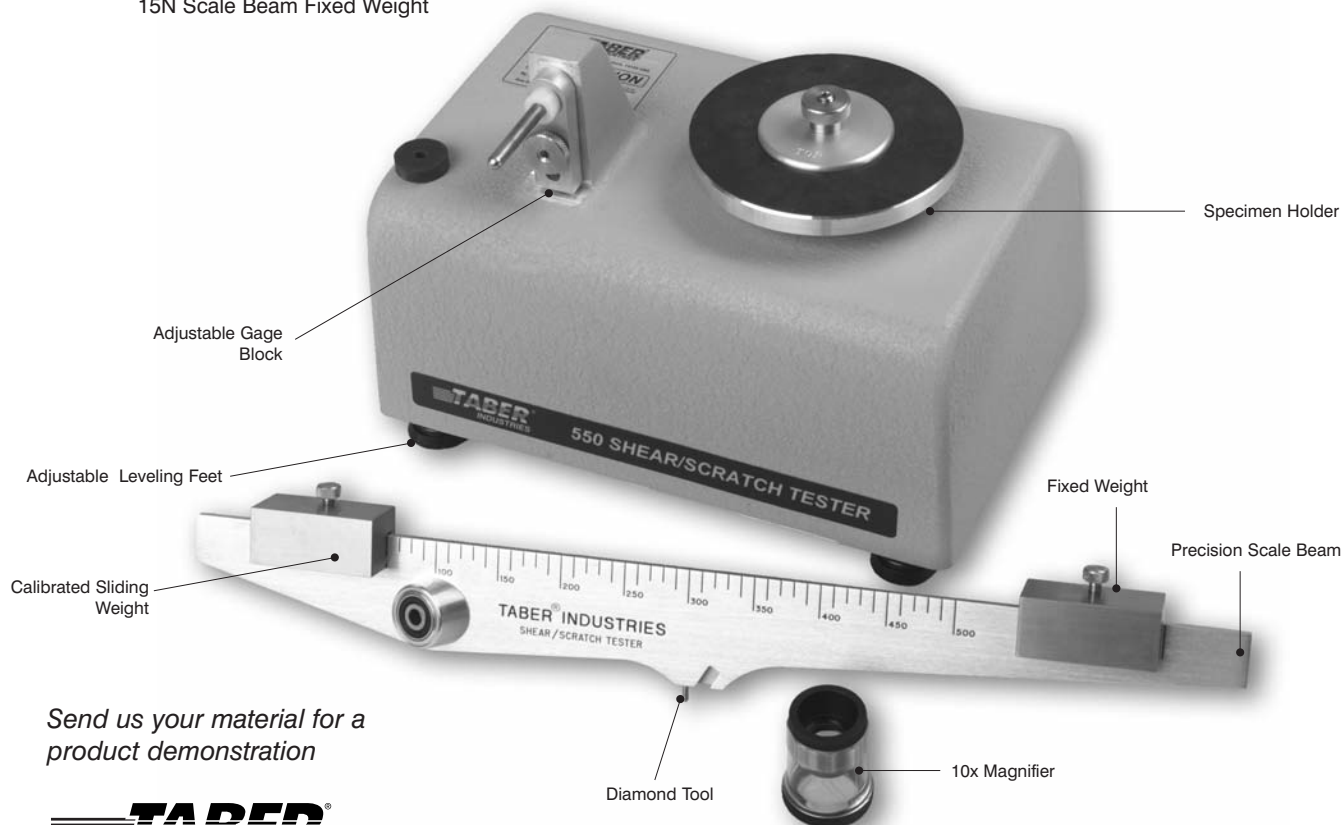

*Send us your material for a
product demonstration*

TABER
INDUSTRIES

Setting the Industry Standard

455 Bryant Street
North Tonawanda, New York 14120 USA

For more information:
Call 1.800.333.5300 (U.S.) or 716.694.4000 (outside U.S.)
Fax 716.694.1450
Email us at sales@taberindustries.com
Visit us at www.taberindustries.com