

# Quiet Cabinet

**TABER**<sup>®</sup>  
INDUSTRIES  
ISO 9001:2000 CERTIFIED

## Taber<sup>®</sup> Abraser Accessory

### **Reduce noise levels associated with the Taber<sup>®</sup> Rotary Abraser operation.**

The typical sound level of operating a Taber Abraser and its vacuum system can reach 78dB. Utilizing sound deadening materials, the Quiet Cabinet reduces this level up to 20%. The result is a noise level comparable to a standard office environment, and enables a normal conversation to be conducted in the vicinity of the testing.

### **Designed for optimal performance.**


The upper cabinet features a hinged Plexiglas viewing window to monitor testing. A detachable front permits easy transfer of the Abraser into or out of the cabinet. Plus, side vacuum hose openings (left or right) allow the top to be used as a standalone cabinet. The lower cabinet holds the Abraser Vacuum unit and includes a built-in exhaust system to provide the proper air circulation. Ample room is available to store test specimens, supplies, and accessories.

### **Built for laboratory or production settings.**

The Quiet Cabinet provides a convenient, dust-free workspace for your Taber Abraser. Constructed with solid wood, a high-pressure laminate outer shell is used for added durability. Aesthetically pleasing, the Quiet Cabinet blends with any laboratory or manufacturing environment.

### **Special top enclosure for Grit Feeder Attachment.**

To accommodate the grit feeder, a large top enclosure is available. Intended for single head Taber Rotary Abrasers and either model 155 or 255 Grit Feeder, the cabinet is 82.5 cm (32.5") in height and features two hinged viewing doors.

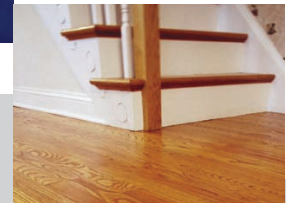

# Quiet Cabinet

Taber® Abraser Accessory

## Available in Two Models:

120V, 60Hz exhaust system  
230V, 50Hz exhaust system

CE

## Standard Features:

Hinged Plexiglas viewing window  
Lower cabinet exhaust system  
Detachable front  
Solid wood construction  
Side vacuum hose openings  
Adjustable height leveling feet  
High-pressure laminate outer shell  
Magnetic door latch


## Quiet Cabinet Specifications

Complete: 76 cm wide x 70 cm deep x 135 cm high  
[30" wide x 24" deep x 53" high]

Top: 76 cm wide x 70 cm deep x 50 cm high  
[30" wide x 24" deep x 19.5" high]

Bottom: 76 cm wide x 70 cm deep x 85 cm high  
[30" wide x 24" deep x 33.5" high]

NOTE: Available as complete unit or Top / Bottom Cabinets may be purchased separately.

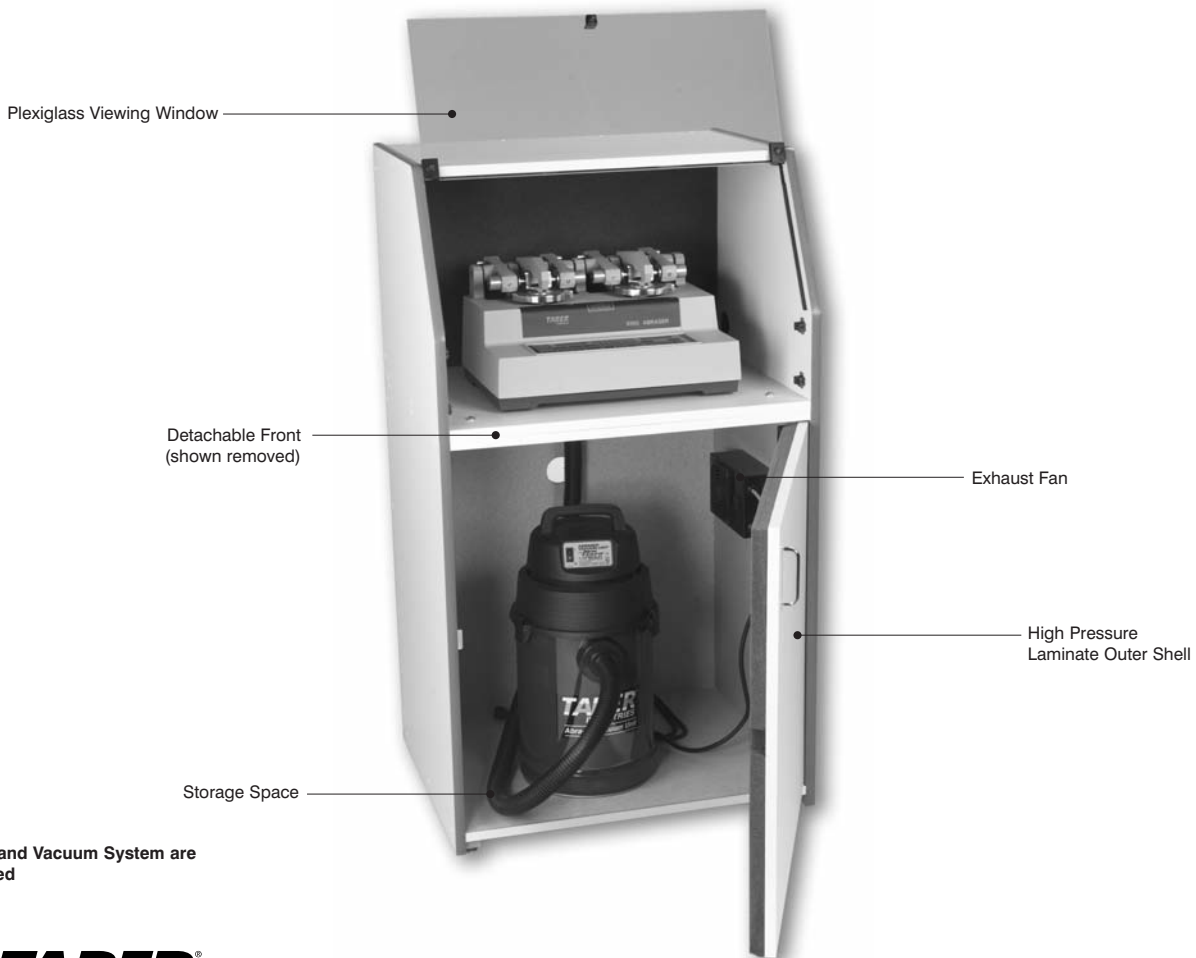

\* Abraser and Vacuum System are not included

**TABER**  
INDUSTRIES

Setting the Industry Standard

455 Bryant Street  
North Tonawanda, New York 14120 USA

For more information:  
Call 1.800.333.5300 (U.S.) or 716.694.4000 (outside U.S.)  
Fax 716.694.1450  
Email us at [sales@taberindustries.com](mailto:sales@taberindustries.com)  
Visit us at [www.taberindustries.com](http://www.taberindustries.com)